

107

-14

Acetaminophen

Salicylamide

Excedrin TabletsColumn: μ Bondapak C₁₈

Solvent: 10 : 90

MeOH : 1% AcOH
acetic acid


Flow rate: 2.0 ml/min

Detector: U.V., 254nm, 0.16 μ ves

1 Tablet in 100 ml MeOH

Caffeine

Aspirin

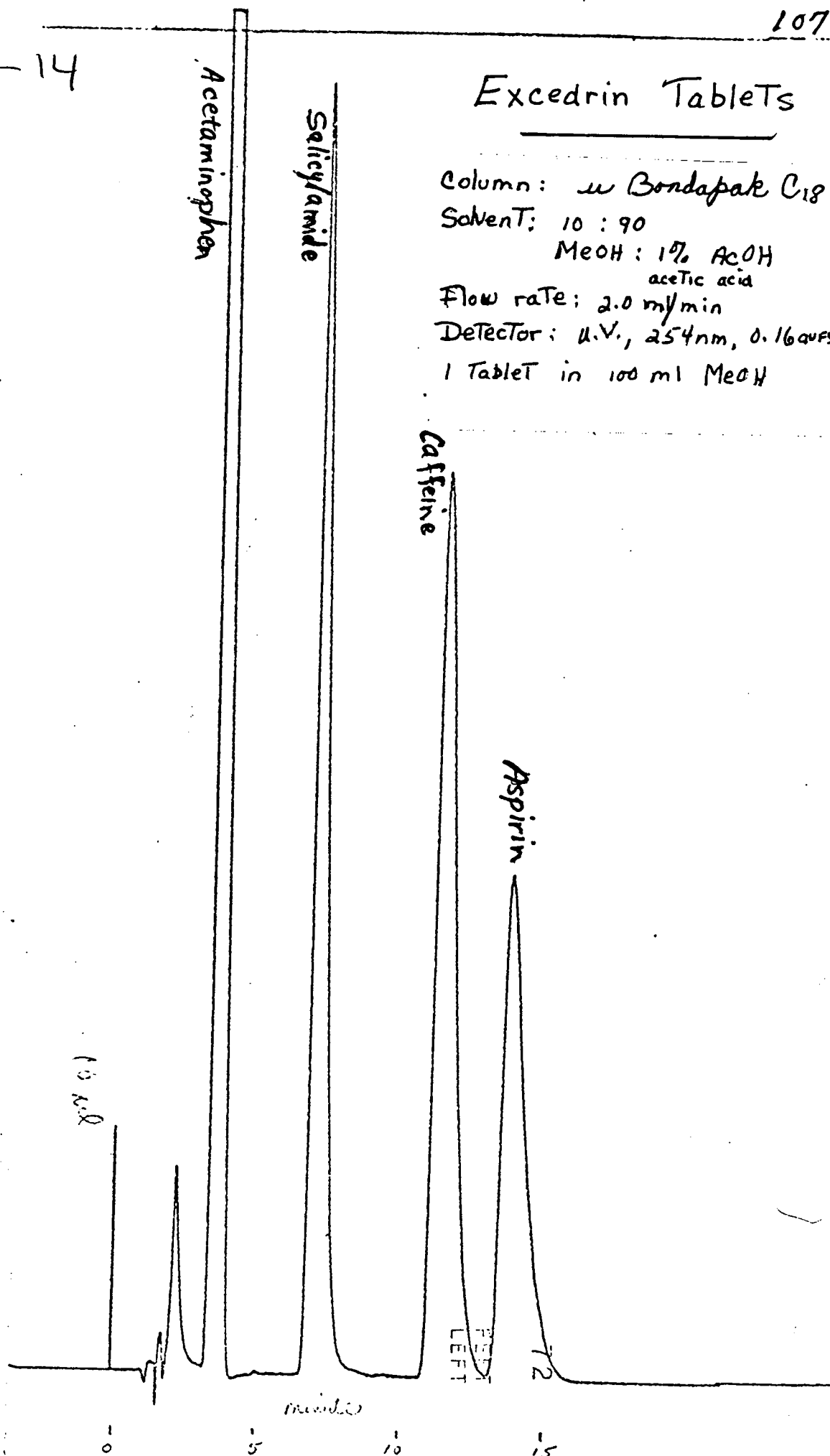

Summary

- **Provides Qualitative and Quantitative Results**
 - **Spectral data for peak identification and peak homogeneity information**
 - **Quantitation of major and minor peaks**
- **Sensitivity Equivalent to Variable Wavelength UV/visible Detectors**
 - **Quantitation of peak less than 0.1%**
 - **Detection of low concentration of coeluting impurities**
- **Compound Confirmation Assures Quality Data**
 - **Peak Purity confirms no coelution detected**
 - **Library Match confirms sample spectra match reference compound spectra**