



FUNGICIDES IN ORANGE JUICE

HPLC METHOD


COLUMN: XTerra® RP₁₈, 4.6 x 100 mm, 3.5 µm
PART NUMBER: 186000438
MOBILE PHASE: 72.5% phosphate buffer, (20 mM pH 6.8)
 27.5% acetonitrile
FLOW RATE: 1.0 mL/min
INJECTION VOLUME: 20 µL
DETECTION: PDA (288nm)

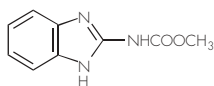

CARBENDAZIM

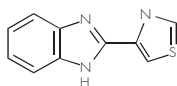


THIABENDAZOLE

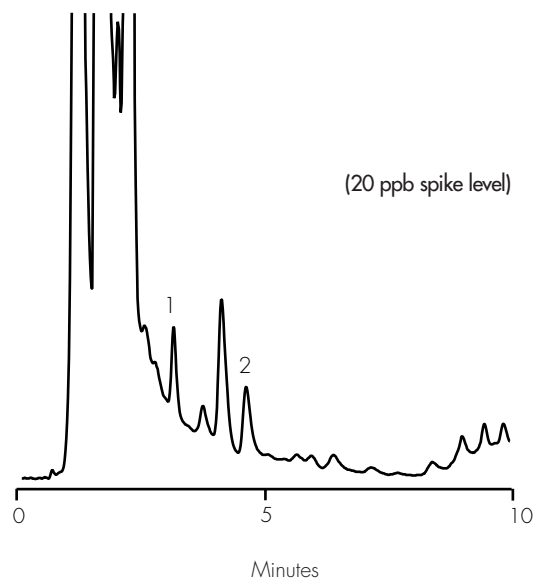

OASIS® MCX SPE METHOD

Conditions for Oasis® MCX Cartridge, 6 cc/150 mg
 Part Number 186000256

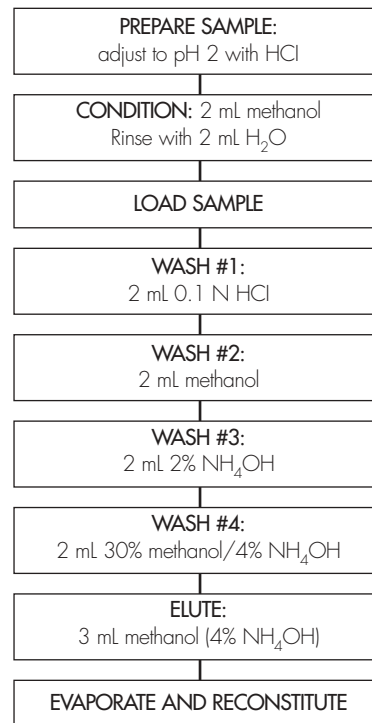

COMPOUNDS	% Recovery, n=4
1. Carbendazim	91% (2% RSD)
2. Thiabendazole	94% (3% RSD)