


Verapamil and Norverapamil

HPLC Method


Column: SymmetryShield™ RP₈, 3.9 x 150 mm, 5 μm
Part number: WAT200655
Mobile phase: 50 mM phosphate pH 7.0/acetonitrile/methanol 41:37:22
Flow rate: 1.0 mL/min
Injection volume: 40 μL
Detection: UV @ 230 nm

USP Tailing Factors


1. 1.10
2. 1.06
3. 1.07

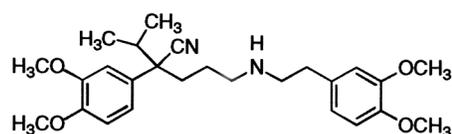

1. Norverapamil

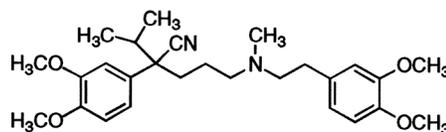

2. Verapamil

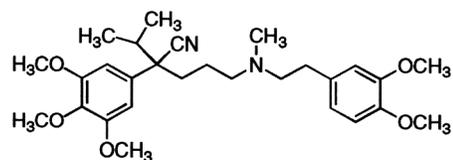

3. Methoxyverapamil (I.S.)

